

# FITTING GUIDELINES

Follow these easy steps to fit successfully

## 1. SELECT INITIAL CONTACT LENS

Upon completion of the refraction, use the **Spectacle to Contact Lens Power Grid\*\*** chart on the back to select your initial trial contact lens power.

### Contact Lens Axis:

- Select the axis for DAILIES® AquaComfort Plus® Toric contact lenses closest to the spectacle axis – 10°, 20°, 70°, 80°, 90°, 100°, 110°, 160°, 170°, 180°.
- Allow the contact lens to settle on the eye.

## 2. EVALUATE CONTACT LENS FIT

### a) Assess centration:

A well-fitted contact lens will have full corneal coverage with good centration.


### b) Assess contact lens movement:

A well-fitted contact lens will show sufficient movement to allow tear exchange under the contact lens during blink in primary or upgaze.

### c) Assess push-up test:

Using lower eyelid, push the contact lens upward. A well-fitted contact lens is easy to displace and returns smoothly to the centered position when released.

## 3. EVALUATE INITIAL CONTACT LENS ORIENTATION

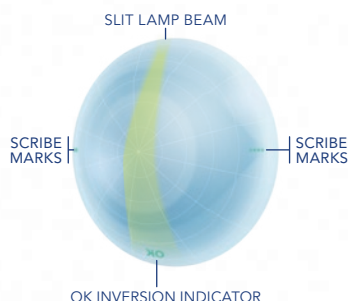

### Tips for viewing the scribe lines:

- Narrow the slit lamp beam.
- Confirm the light beam passes through the center of the pupil. Scribe lines are located at 3 and 9 o'clock. When fitting, the straight line of the "K" in the OK Inversion Indicator can also be used as the 6 o'clock position to measure rotation.
- Determine the amount of rotation, if any, and compensate accordingly.

**NOTE:** OK Inversion Indicator may orient in either the 6 or 12 o'clock position

## 4. EVALUATE VISUAL ACUITY

Check visual acuity and perform a spherical over-refraction if necessary. If the spherical over-refraction does not yield acceptable vision, perform a spherocylindrical over-refraction. For the correct contact lens power from this over-refraction, **call consultation at 1-800-241-7468 for additional support.** Visit [myalcon.com](http://myalcon.com) for more details.



PERFORMANCE DRIVEN BY SCIENCE™

# SPECTACLE TO CONTACT LENS POWER GRID

Vertex Correction: 12mm

## SPECTACLE CYLINDER POWERS

SPECTACLE SPHERE POWERS

	-0.75		-1.00		-1.25		-1.50		-1.75		-2.00		-2.25		-2.50	
	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl
-9.00	-8.50	sph	-8.00	-0.75	-8.00	-0.75	-8.00	-1.25	-8.00	-1.25	-8.00	-1.75	-8.00	-1.75	-8.00	-1.75
-8.75	-8.00	sph	-8.00	-0.75	-8.00	-0.75	-8.00	-1.25	-8.00	-1.25	-8.00	-1.75	-8.00	-1.75	-8.00	-1.75
-8.50	-8.00	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-1.75
-8.25	-7.75	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-1.75
-8.00	-7.75	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-1.75
-7.75	-7.50	sph	-7.00	-0.75	-7.00	-0.75	-7.00	-1.25	-7.00	-1.25	-7.00	-1.75	-7.00	-1.75	-7.00	-1.75
-7.50	-7.25	sph	-7.00	-0.75	-7.00	-0.75	-7.00	-1.25	-7.00	-1.25	-7.00	-1.75	-7.00	-1.75	-7.00	-1.75
-7.25	-7.00	sph	-6.50	-0.75	-6.50	-0.75	-6.50	-1.25	-6.50	-1.25	-6.50	-1.75	-6.50	-1.75	-6.50	-1.75
-7.00	-6.75	sph	-6.50	-0.75	-6.50	-0.75	-6.50	-1.25	-6.50	-1.25	-6.50	-1.75	-6.50	-1.75	-6.50	-1.75
-6.75	-6.50	sph	-6.00	-0.75	-6.00	-0.75	-6.00	-1.25	-6.00	-1.25	-6.00	-1.75	-6.00	-1.75	-6.00	-1.75
-6.50	-6.00	-0.75	-6.00	-0.75	-6.00	-0.75	-6.00	-1.25	-6.00	-1.25	-6.00	-1.75	-6.00	-1.75	-6.00	-1.75
-6.25	-5.75	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-1.75
-6.00	-5.50	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-1.75
-5.75	-5.50	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-1.75
-5.50	-5.25	-0.75	-5.25	-0.75	-5.25	-0.75	-5.25	-1.25	-5.25	-1.25	-5.25	-1.75	-5.25	-1.75	-5.25	-1.75
-5.25	-5.00	-0.75	-5.00	-0.75	-5.00	-0.75	-5.00	-1.25	-5.00	-1.25	-5.00	-1.75	-5.00	-1.75	-5.00	-1.75
-5.00	-4.75	-0.75	-4.75	-0.75	-4.75	-0.75	-4.75	-1.25	-4.75	-1.25	-4.75	-1.75	-4.75	-1.75	-5.00	-1.75
-4.75	-4.50	-0.75	-4.50	-0.75	-4.50	-0.75	-4.50	-1.25	-4.50	-1.25	-4.50	-1.75	-4.50	-1.75	-4.75	-1.75
-4.50	-4.25	-0.75	-4.25	-0.75	-4.50	-0.75	-4.25	-1.25	-4.25	-1.25	-4.25	-1.75	-4.25	-1.75	-4.50	-1.75
-4.25	-4.00	-0.75	-4.00	-0.75	-4.25	-0.75	-4.00	-1.25	-4.00	-1.25	-4.00	-1.75	-4.00	-1.75	-4.25	-1.75
-4.00	-4.00	-0.75	-4.00	-0.75	-4.00	-0.75	-4.00	-1.25	-4.00	-1.25	-4.00	-1.75	-4.00	-1.75	-4.25	-1.75
-3.75	-3.75	-0.75	-3.75	-0.75	-3.75	-0.75	-3.75	-1.25	-3.75	-1.25	-3.75	-1.75	-3.75	-1.75	-3.75	-1.75
-3.50	-3.50	-0.75	-3.50	-0.75	-3.50	-0.75	-3.50	-1.25	-3.50	-1.25	-3.50	-1.75	-3.50	-1.75	-3.50	-1.75
-3.25	-3.25	-0.75	-3.25	-0.75	-3.25	-0.75	-3.25	-1.25	-3.25	-1.25	-3.25	-1.75	-3.25	-1.75	-3.25	-1.75
-3.00	-3.00	-0.75	-3.00	-0.75	-3.00	-0.75	-3.00	-1.25	-3.00	-1.25	-3.00	-1.75	-3.00	-1.75	-3.00	-1.75
-2.75	-2.75	-0.75	-2.75	-0.75	-2.75	-1.25	-2.75	-1.25	-2.75	-1.25	-2.75	-1.75	-2.75	-1.75	-2.75	-1.75
-2.50	-2.50	-0.75	-2.50	-0.75	-2.50	-1.25	-2.50	-1.25	-2.50	-1.25	-2.50	-1.75	-2.50	-1.75	-2.75	-1.75
-2.25	-2.25	-0.75	-2.25	-0.75	-2.25	-1.25	-2.25	-1.25	-2.25	-1.75	-2.25	-1.75	-2.25	-1.75	-2.50	-1.75
-2.00	-2.00	-0.75	-2.00	-0.75	-2.00	-1.25	-2.00	-1.25	-2.00	-1.75	-2.00	-1.75	-2.00	-1.75	-2.25	-1.75
-1.75	-1.75	-0.75	-1.75	-0.75	-1.75	-1.25	-1.75	-1.25	-1.75	-1.75	-1.75	-1.75	-2.00	-1.75	-2.00	-1.75
-1.50	-1.50	-0.75	-1.50	-0.75	-1.50	-1.25	-1.50	-1.25	-1.50	-1.75	-1.50	-1.75	-1.75	-1.75		
-1.25	-1.25	-0.75	-1.25	-0.75	-1.25	-1.25	-1.25	-1.25	-1.25	-1.75	-1.25	-1.75	-1.50	-1.75		
-1.00	-1.00	-0.75	-1.00	-0.75	-1.00	-1.25	-1.00	-1.25	-1.00	-1.75	-1.00	-1.75	-1.25	-1.75		
-0.75	-0.75	-0.75	-0.75	-0.75	-0.75	-1.25	-0.75	-1.25	-0.75	-1.75	-0.75	-1.75	-1.00	-1.75		
-0.50	-0.50	-0.75	-0.50	-0.75	-0.50	-1.25	-0.50	-1.25	-0.50	-1.75	-0.50	-1.75	-0.75	-1.75		
-0.25	-0.25	-0.75	-0.25	-0.75	-0.25	-1.25	-0.25	-1.25	-0.25	-1.75	-0.25	-1.75	-0.50	-1.75		
+0.00	+0.00	-0.75	+0.00	-0.75	+0.00	-1.25	+0.00	-1.25	+0.00	-1.75	+0.00	-1.75	-0.25	-1.75		
+0.25	+0.25	-0.75	+0.25	-0.75	+0.25	-1.25	+0.25	-1.25	+0.25	-1.75	+0.25	-1.75	0.00	-1.75		
+0.50	+0.50	-0.75	+0.50	-0.75	+0.50	-1.25	+0.50	-1.25	+0.50	-1.75	+0.50	-1.75	+0.25	-1.75		
+0.75	+0.75	-0.75	+0.75	-0.75	+0.75	-1.25	+0.75	-1.25	+0.75	-1.75	+0.75	-1.75	+0.50	-1.75		
+1.00	+1.00	-0.75	+1.00	-0.75	+1.00	-1.25	+1.00	-1.25	+1.00	-1.75	+1.00	-1.75	+0.75	-1.75		
+1.25	+1.25	-0.75	+1.25	-0.75	+1.25	-1.25	+1.25	-1.25	+1.25	-1.75	+1.25	-1.75	+1.00	-1.75		
+1.50	+1.50	-0.75	+1.50	-0.75	+1.50	-1.25	+1.50	-1.25	+1.50	-1.75	+1.50	-1.75	+1.25	-1.75		
+1.75	+1.75	-0.75	+1.75	-0.75	+1.75	-1.25	+1.75	-1.25	+1.75	-1.75	+1.75	-1.75	+1.50	-1.75		
+2.00	+2.00	-0.75	+2.00	-0.75	+2.00	-1.25	+2.00	-1.25	+2.00	-1.75	+2.00	-1.75	+1.75	-1.75		
+2.25	+2.25	-0.75	+2.25	-0.75	+2.25	-1.25	+2.25	-1.25	+2.25	-1.75	+2.25	-1.75	+1.75	-1.75		
+2.50	+2.50	-0.75	+2.50	-0.75	+2.50	-1.25	+2.50	-1.25	+2.50	-1.75	+2.50	-1.75	+2.25	-1.75		
+2.75	+2.75	-0.75	+2.75	-0.75	+2.75	-1.25	+2.75	-1.25	+2.75	-1.75	+2.75	-1.75	+2.50	-1.75		
+3.00	+3.00	-0.75	+3.00	-0.75	+3.00	-1.25	+3.00	-1.25	+3.00	-1.75	+3.00	-1.75	+2.75	-1.75		
+3.25	+3.25	-0.75	+3.25	-0.75	+3.25	-1.25	+3.25	-1.25	+3.25	-1.75	+3.25	-1.75	+3.00	-1.75		
+3.50	+3.50	-0.75	+3.50	-0.75	+3.50	-1.25	+3.50	-1.25	+3.50	-1.75	+3.50	-1.75	+3.25	-1.75		
+3.75	+3.75	-0.75	+3.75	-0.75	+3.75	-1.25	+3.75	-1.25	+3.75	-1.75	+3.75	-1.75	+3.50	-1.75		
+4.00	+4.00	-0.75	+4.00	-0.75	+4.00	-1.25	+4.00	-1.25	+4.00	-1.75	+4.00	-1.75	+3.75	-1.75		
	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl

TRIAL CONTACT LENS POWERS

HOW TO USE THE **DAILIES® AQUACOMFORT PLUS® TORIC** CONTACT LENS SELECTION GRID  
**NOTE:** There is no need to correct for vertex distance, as the chart will do this for you.

1. Locate the sphere power in the wearer's spectacle prescription on the left side labeled "Spectacle Sphere Powers."
2. Locate the DAILIES® AquaComfort Plus® Toric trial contact lens power underneath the spectacle minus cylinder power in the wearer's spectacle Rx.
3. For trial contact lenses, select the DAILIES® AquaComfort Plus® Toric contact lens axis closest to the refractive or spectacle prescription axis.
4. When ordering the final contact lens, account for contact lens rotation as appropriate.

# REFITTING FROM FOCUS® DAILIES® TORIC CONTACT LENSES GUIDELINES



**FOCUS® DAILIES® TORIC**  
CONTACT LENSES ARE AVAILABLE IN  
THE FOLLOWING PARAMETERS:

<b>BASE CURVE (mm):</b>	8.6
<b>DIAMETER (mm):</b>	14.2
<b>SPHERE POWERS (D):</b>	+4.00 to -6.00 (0.25D steps) -6.50 to -8.00 (0.50D steps)
<b>CYLINDER POWERS (D):</b>	-0.75; -1.50
<b>AXES:</b>	20°; 70°; 90°; 110°; 160°; 180°



**DAILIES® AQUACOMFORT PLUS® TORIC**  
CONTACT LENSES ARE AVAILABLE IN  
THE FOLLOWING PARAMETERS:

<b>BASE CURVE (mm):</b>	8.8
<b>DIAMETER (mm):</b>	14.4
<b>SPHERE POWERS (D):</b>	+4.00 to -6.00 (0.25D steps) -6.50 to -8.00 (0.50D steps)
<b>CYLINDER POWERS (D):</b>	-0.75; -1.25; -1.75
<b>AXES:</b>	10°; 20°; 70°; 80°; 90°; 100°; 110°; 160°; 170°; 180°

## IT'S EASY TO MAKE THE SWITCH

### EXPANDED AXIS OPTIONS

- All current FOCUS® DAILIES® Toric cylinder axes are also available in DAILIES® AquaComfort Plus® Toric contact lenses.
  - If a wearer's cylinder axis is already optimally corrected with their prescribed FOCUS® DAILIES® Toric lens—based on spectacle Rx and any observed rotation—it's not necessary to change the cylinder axis when refitting into DAILIES® AquaComfort Plus® Toric lenses.
  - If one of the additional DAILIES® AquaComfort Plus® Toric lens axes would more optimally correct a wearer's astigmatic refractive error—based on spectacle Rx or any observed lens rotation—you may opt to trial-fit in one of these axes.

### EXPANDED CYLINDER POWER OPTIONS

- For current FOCUS® DAILIES® Toric wearers in the -1.50D cylinder power, determine if the -1.25D or -1.75D cylinder power is more appropriate.
  - Select the DAILIES® AquaComfort Plus® Toric cylinder power that is closer to the vertex-corrected cylinder power.

For additional information or assistance with refitting FOCUS® DAILIES® Toric contact lens wearers into DAILIES® AquaComfort Plus® Toric contact lenses, please contact Alcon technical consultation at **1-800-241-7468**.

## DAILIES® AquaComfort Plus® Toric contact lenses product information

<b>MATERIAL:</b>	nelfilcon A
<b>MOISTURIZING AGENTS:</b>	Unique Blink-Activated Moisture System <ul style="list-style-type: none"> <li>• Hydroxypropyl methylcellulose (HPMC)</li> <li>• Polyethylene glycol (PEG)</li> <li>• Polyvinyl alcohol (PVA)</li> </ul>
<b>WATER CONTENT:</b>	69%
<b>HANDLING TINT:</b>	Light blue VISITINT®
<b>SPHERE POWERS (D):</b>	+4.00 to -6.00 (0.25D steps) -6.50 to -8.00 (0.50D steps)
<b>CYLINDER POWERS (D):</b>	-0.75; -1.25; -1.75
<b>AXES:</b>	10°; 20°; 70°; 80°; 90°; 100°; 110°; 160°; 170°; 180°
<b>DIAMETER (mm):</b>	14.4
<b>BASE CURVE (mm):</b>	8.8
<b>CENTER THICKNESS (mm):</b>	0.10 @ -3.00D
<b>SCRIBE MARKS:</b>	Scribe marks at 3 and 9 o'clock
<b>INVERSION INDICATOR:</b>	OK
<b>DESIGN:</b>	Precision Curve™ Lens Design/Toric back surface design
<b>WEAR/REPLACEMENT SCHEDULE:</b>	Daily Wear, Daily Disposable
<b>PACK SIZES:</b>	Trials – 5 lenses / Revenue – 30 and 90 lens packs

For further questions, or if the visual acuity is still not satisfactory, please call consultation at **1-800-241-7468** for additional support.

Visit **myalcon.com** for more details.



\*Results may vary. See USA package insert for details. \*\*Determination of which final cylinder axis to order will be made after compensation for contact lens orientation.

See product instructions for complete wear, care, and safety information.

© 2013 Novartis 8/13 DAT13018MS



**Alcon**  
a Novartis company